



SafeLock™

Automatic needle guard for pre-filled syringe

Device Type

Disposable Safety Syringe

Primary Container

1mL Long Syringe with Staked Needle or 2.25mL Syringe with Staked Needle

Dosage Type

Fixed | Full Contents

Fill Volume

0.1-1mL | 1.1-2.25mL

Route of Administration

Subcutaneous

Handling

Manual, one-handed use, automatic lock-out at end of dose

Injection Feedback

Audible click at beginning & end of injection, viewing window


Handling Instructions



1
Remove cap



2
Press & hold against injection site vertically

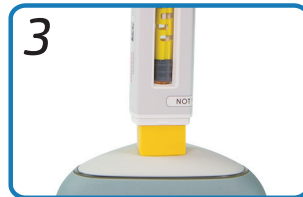

3
Press down plunger rod to administer the drug, causing needle cover to lock into place



4
Discard in sharps container

Key Features

- ✓ Simple 3-step operation
- ✓ Audible & visual feedback
- ✓ One-handed operation
- ✓ Compact design
- ✓ Locking needle cover
- ✓ Large labeling area
- ✓ Preconfigured, with colors customized to client specifications
- ✓ Compatible with Rigid Needle Shield